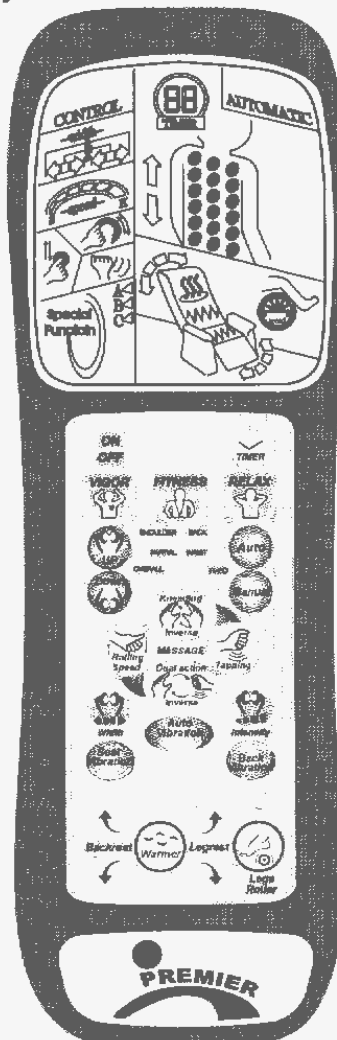



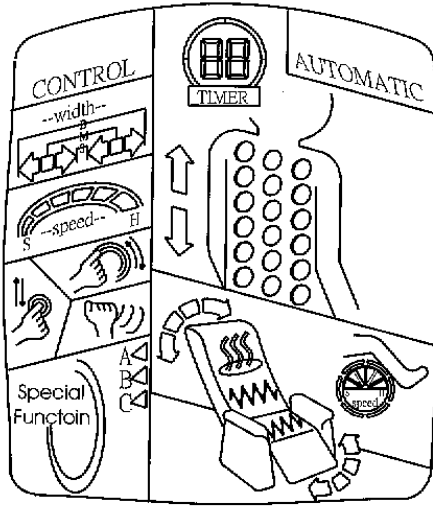




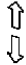











Product Function And Descriptions

※Please read the instruction manual carefully in order to operate the chair proficiently. Refer to below for the function keys.



- DO NOT use foreign or hard objects to press the controller's buttons, including nails. Improper use will result in damaging the controller.
- DO NOT drop the controller. Keep hands dry to use the controller.
- DO NOT sit or step on the controller.
- DO NOT press the buttons too hard.
- TURN OFF the power AFTER finishing using the chair.

LCD Features

 Automatic	 <p>LCD Display Screen</p>	
 Timer		
 Positions of massage rollers		
 Rollers' width adjustment		
 Speed adjustment		
 Direction of massage rollers	 Rolling	
 Back warmer	 Kneading/Tapping	
 Back vibration massage	 Special function	
 Seat vibration massage	 Foot rolling	
 Kneading	 Backrest Recline	
 Tapping	 Legrest position	



Power (Start / Stop)

Press the button to turn on/off the chair. The LCD shows "00" when the chair is on. When using the chair, press the button to reset the chair to fully upright position, and return to standby "00" mode.



Timer

Press the button to set the massage time from 5 to 30 minutes with 5-minute interval.



Special Function: VIGOR

Press the button to turn on/off the 30-minute programmed massage, and the LCD shows the resting time and the functioning massage.

Before massage, the massage chair would turn on automatic shoulder detection. This is to locate the acupuncture points of the user and provide the most comfortable massage. The function would automatically adjust the width and speed of massage rollers. When the roller is upper the shoulder, the manual adjustment for width is not available. (See LCD features)



Special Function: FITNESS

Press the button to turn on/off the function. The 15-minute program is preset. The LCD would show the progress and the time left. The programmed massage would automatically adjust the width and the speed of roller for the user.

Before massage, the massage chair would turn on automatic shoulder detection. This is to locate the acupuncture points of the user and provide the most comfortable massage.

The function would automatically adjust the width and speed of massage rollers. When the roller is upper the shoulder, the manual adjustment for width is not available. (See LCD features)



Special Function: RELAX

Press the *button* to turn on/of the function. The preset time is 10 minutes, but it is manually adjustable up to 99 minutes. The LCD would show the progress and the time left. In the first 3-minute *massage*, the *programmed* *massage* would automatically adjust the width and the speed of roller for the user. After the first 3-minute *massage*, the *massage* would be off automatically, so the user can have a short rest. In the last 3 minutes, the *massage chair* would function to wake up the user. (See LCD features)



Roller Positioning Upward

When using the manual function, press and hold the *button* to adjust positions of the roller upward to the desired position, and then, release the *button*. (Not available for Special and auto functions)



Roller Positioning Downward

When using the manual function, press and hold the *button* to adjust positions of the rollers downward to the desired position, and then, release the *button*. (Not available for Special and Auto functions)



Auto

Press the *button* to turn on the function for whole, partial back, shoulders, waist, or fixed point. The preset time is 15 minutes, shown on the LCD screen. The flashing LED indicates the working function user selected. The massaging width and speed are adjusted automatically by the preset program. Before *massage*, the *massage chair* would turn on automatic shoulder detection. This is to locate the acupuncture points of the user and provide the most comfortable *massage*. The function would automatically adjust the width and speed of *massage* rollers. When the roller is upper the shoulder, the manual adjustment for width is not available. (See LCD features)



MANUAL

Press the button to select the functions for whole back, partial, shoulders, waist, to fixed point. The preset time is 15 minutes, and it is adjustable up to 30 minutes. The selected functions are also shown on the panel through LED flashing. The user can also select different function at once.




Kneading

Press the button to activate Kneading/Inverse Kneading/off the function.

One click the button leads the kneading and inward rotating massage and duple-click does the kneading and ourward rotating massage. User can select the different functions at once. The selected function would be shown by the flashing LED. While the kneading is selected, user can adjust the speed of roller by press the ROLLING/SPEED. There are 6 speeds available. The selected speed is shown on the LCD with 3-second display. (See LCD features)



Rolling/Speed

Press the button to turn on/off the chiropractice-like massage function with the time setting from 5 to 30 minutes, and adjust the up-and-down speed of rollers (not available for fixed massage). Pressing "  " could

adjust the speed as well, and the speed would be shown on the LCD with 3-second display. The rolling width would be narrowed automatically as the roller moves upper the shoulder if the WIDE rolling is set. The rolling width remains the same if the MEDIUM or NARROW is set.



Tapping

Press the button to turn on/off the function. User can select different manual massages at once, and set the time from 5 to 30 minutes. While the tapping is selected, user can adjust the speed of roller by press the ROLLING/SPEED. There are 6 speeds available. The selected speed is shown on the LCD with 3-second display. The manual adjustment to narrow the roller is not available when the roller is upper the shoulder. (See LCD features)



Dual Action

Press the button to activate Dual Action/Inverse Dual Action /off the function combining the tapping and kneading massage, and the time is adjustable from 5 to 30 minutes.

Clicking the button once leads the kneading and inward rotating massage and double-click does the kneading and outward rotating massage.

While the kneading is selected, user can adjust the speed of roller by pressing the ROLLING/SPEED. There are 6 speeds available. The selected speed is shown on the LCD with 3-second display. (See LCD features)



Width Adjustment

Width adjustment for tapping and rolling. User can select wide, medium or narrow for personal comfort ability. It is not available for the kneading and dual functions, and the Special and Auto functions. The manual adjustment for narrowing is not available when the roller is upper the shoulder.

※Not functionable when Auto Programs are activated.



Speed Adjustment

There are 6 speeds available for tapping, kneading, and dual action. It is also available for rolling massage by clicking the Roller Positioning Upward/Downward or the ROLLING/SPEED. When the Special or Auto function is selected, the speed is controlled by the default massage programs.

※Not functionable when Auto Programs are activated.



Seat Vibration

Press the button to turn on/off the function. Press again to select the strong/weak vibrating intensity. (See the LCD features. The function is not available for Special and Auto functions.)



Auto Vibration

Press the button to turn on/off the function. Press again to select the strong/weak vibrating intensity. (See the LCD features. The function is not available for Special and Auto functions.)



Back Vibration

Press the button to turn on/off the function. Press again to select the strong/weak vibrating intensity. (See the LCD features. The function is not available for Special and Auto functions.)



Backrest Adjustment

Press and hold the up or down arrow to adjust the angle of the backrest. Double-click can synchronize the movement of the backrest and footrest. (The function is available even though the power of remote controller is off.)



Back Warmer

Press the button to turn on/off the warming device in the backrest. The device can work with other functions. (Please refer to the Instruction on Page 04)



Footrest Adjustment

Press the button to control the footrest upward and downward. Double-click can synchronize the movement of backrest and footrest.



Legs Roller

Press the button to turn on/off the function, or select the weak / medium / strong intensity to massage. (See the display for intensit on the LCD)

Product Functions And Descriptions

Instruction on Using the Manual Calf Air Massage

- Find the hand air pump
- Squeeze the air pump to inflate the calf airbags.
- Press button B to deflate calf airbags, and place the hand air pump inside the pocket at the rear of the left armrest when not in use.

